

TRIPLEBODY

ADVANCED PROFESSIONAL SYSTEM



"TRIPLEBODY" is a 3 in 1 innovative beauty device that has MF (Middle Frequency), RF and Laser functions in 1 device.



Creative Approach !

Innovative Design !



- ☐ Operates MF + RF + Laser Simultaneously!
- ☐ 8 Selective Pads that you can apply Anywhere You Want!
- ☐ Simultaneous Treatment with 8 pads!
- ☐ Individual Waveforms and Level Adjustments on each 8 pads!
- ☐ Build-in Temperature Sensor – Automatic OFF function!
- ☐ No Side Effect! Easy Treatment!

EFFECTS


MF	RF
<ul style="list-style-type: none"><input type="checkbox"/> Dissolve Fat Cells<input type="checkbox"/> Decrease Cellulites<input type="checkbox"/> Relaxing Muscles<input type="checkbox"/> Detox<input type="checkbox"/> Blood Circulation<input type="checkbox"/> Accelerate Metabolism	<ul style="list-style-type: none"><input type="checkbox"/> Improves Wrinkles<input type="checkbox"/> Anti-Aging<input type="checkbox"/> Skin Elasticity
	Laser
	<ul style="list-style-type: none"><input type="checkbox"/> Dissolve Fat Cells<input type="checkbox"/> Decrease Cellulites



Middle Frequency



6 types of waveforms - User Selective

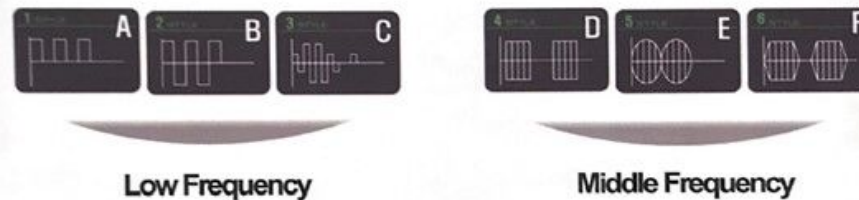

- ☐ 6 types of different waveforms
- ☐ Operation Time : 0~99mins. (default setting 30mins, Max. 99mins.)
- ☐ Level : 0~99
- ☐ Operation time and level are easily adjustable depend on user's convenience.

6 types of waveforms - User Selective



01 |

Lipolysis Current

Make the sharp body shape, low frequency that directly stimulates the part where body type decided.



02 |

Isometric Current

Low frequency current that has special usage on shock range, incident time, modulation.



03 |

Isotonic Current

An alternating current with short impact that almost same with passive motion treatment or recovering muscle.



04 |

Bipolar Interferential Current [Pausing Time Exist]

Produce middle frequency (4KHz). Increases amplitude rhythmically and pause and resume the consecutive operation time, modulates the standard bandwidth of 4KHz.



05 |

Bipolar Interferential Current

Produce middle frequency (4KHz). Increases amplitude rhythmically and modulates the standard bandwidth of 4KHz.



06 |

Russian Current


Russian current was designed by Russia for the athletes and modified from sustained alternating current that feels very relaxed and comfortable.

4 kinds of program - Selective Program

- ☐ Selective and/or combined program among 6 types of waveform.
- ☐ Can select a program out of 4 programs.
- ☐ Default operation times are set by each program.
- ☐ Levels of each program are adjustable within permissible range.

Program	Operation Time (min.)	Level
1	42 mins.	1 ~ 88
2	39 mins.	1 ~ 74
3	40 mins.	1 ~ 62
4	47 mins.	1 ~ 88

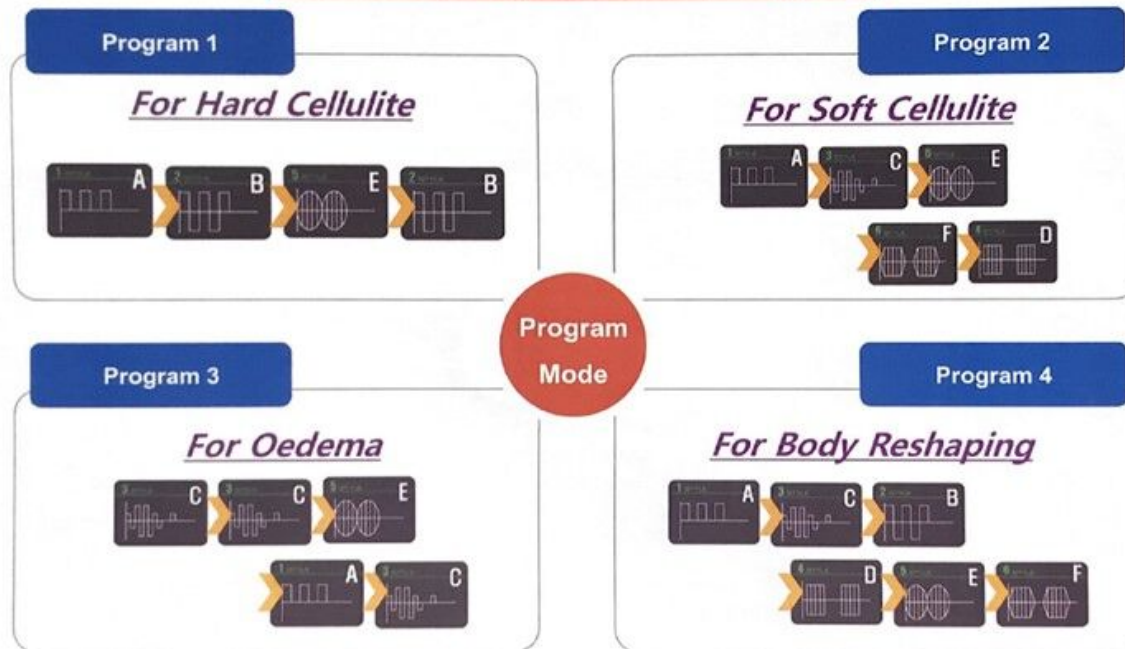

4 kinds of program - **Selective Program**



RF

Transfer the energy into the deep skin through RF to accelerate the collagen creation to improve wrinkle and anti-aging and at the same time, it accelerate the elastin to improve skin elasticity.

- ☐ Operating Frequency : 1 MHz
- ☐ Middle frequency outputs simultaneously on 8 pads
- ☐ Levels of RF are adjustable after operation
- ☐ Level : 5 levels (1 ~ 5)
- ☐ Built-in temperature sensor keeps stable set temperature (Automatic OFF)

Level	Automatic OFF
1	37°C
2	38°C
3	39°C
4	40°C
5	40.5°C



Laser

Laser

When irradiate low output laser on fat area using 635nm diode, the structure of fat cells changes and becomes easy for lipolysis, then it is effective for decreasing fat cells and cellulites.

- ☐ Wave Length : 635nm
- ☐ 5mW x 5 Laser Diodes / PAD
- ☐ Once START, MF & RF output simultaneously
- ☐ Laser Operation : By green light indicator on the pad

Thank you